

交換作業は下記手順で行ってください。

1. 伸展ストッパーゴム交換要領
 - 1) 膝を曲げ、マイナスドライバー等でストッパーゴムの引っかけて抜いてください。(図1参照)
 - 2) 新しいストッパーゴスをマイナスドライバー等ではめ込んでください。
2. 屈曲ストッパーゴム交換要領
 - 1) 膝を完全に伸展させ、マイナスドライバー等で屈曲ストッパーゴムの引っかけて抜いてください。(図2参照)
 - 2) 新しいストッパーゴスをマイナスドライバー等ではめ込んでください。

【Stopper rubber replacement procedure】

Please refer to the following instructions for replacing the extension/flexion stopper rubbers.

1. Extension stopper rubber (See Fig.1.)
 - 1) With the knee bent, using a flat screwdriver or the like, pick up and take out the stopper rubbers.
 - 2) Set new stopper rubbers in place. Make sure to set them completely by using a flat screwdriver or the like.
2. Flexion stopper rubber (See Fig.2.)
 - 1) With the knee fully extended, using a flat screwdriver or the like, pick up and take out the stopper rubbers.
 - 2) Set new stopper rubbers and make sure to set them completely by using a minus driver or the like.

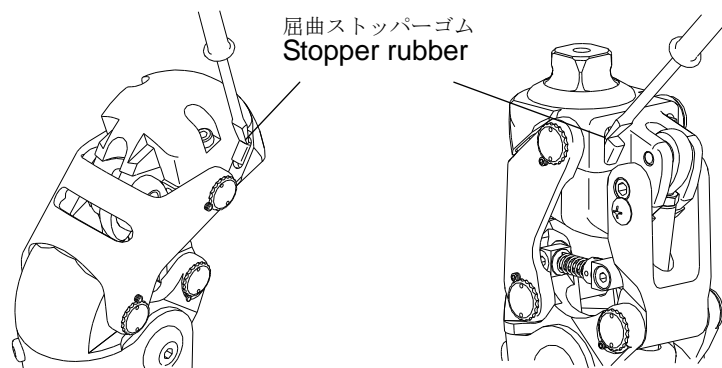


Fig. 1

Fig. 2

お問い合わせ先 ナブテスコ株式会社 住環境カンパニー 福祉事業推進部 営業グループ
 〒658-0024 神戸市東灘区魚崎浜町 35
 TEL:078-413-2724 FAX:078-413-2725

Nabtesco Corporation Accessibility Innovations Company Assistive Products Department
 35, Uozakihama-machi, Higashinada-ku, Kobe, 658-0024, Japan
 PHONE : +81-78-413-2724 FAX : +81-78-413-2725

交換作業は下記手順で行ってください。

3. 伸展ストッパーゴム交換要領
 - 1) 膝を曲げ、マイナスドライバー等でストッパーゴムの引っかけて抜いてください。(図1参照)
 - 2) 新しいストッパーゴスをマイナスドライバー等ではめ込んでください。
4. 屈曲ストッパーゴム交換要領
 - 3) 膝を完全に伸展させ、マイナスドライバー等で屈曲ストッパーゴムの引っかけて抜いてください。(図2参照)
 - 4) 新しいストッパーゴスをマイナスドライバー等ではめ込んでください。

【Stopper rubber replacement procedure】

Please refer to the following instructions for replacing the extension/flexion stopper rubbers.

1. Extension stopper rubber (See Fig.1.)
 - 1) With the knee bent, using a flat screwdriver or the like, pick up and take out the stopper rubbers.
 - 2) Set new stopper rubbers in place. Make sure to set them completely by using a flat screwdriver or the like.
2. Flexion stopper rubber (See Fig.2.)
 - 1) With the knee fully extended, using a flat screwdriver or the like, pick up and take out the stopper rubbers.
 - 2) Set new stopper rubbers and make sure to set them completely by using a minus driver or the like.

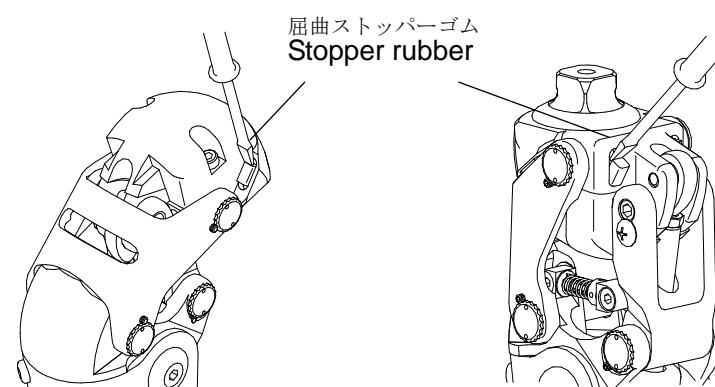


Fig. 1

Fig. 2

お問い合わせ先 ナブテスコ株式会社 住環境カンパニー 福祉事業推進部 営業グループ
 〒658-0024 神戸市東灘区魚崎浜町 35
 TEL:078-413-2724 FAX:078-413-2725

Nabtesco Corporation Accessibility Innovations Company Assistive Products Department
 35, Uozakihama-machi, Higashinada-ku, Kobe, 658-0024, Japan
 PHONE : +81-78-413-2724 FAX : +81-78-413-2725